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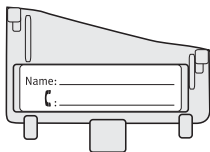
## Declaration of conformity

---

Hereby, Unitron Hearing declares that Smart Control is in compliance with the essential requirements and other relevant provisions of Directive 93/42/EEC. To obtain a copy of the Declaration of Conformity, please contact your local representative from the list of Unitron Hearing Distributors.

This device complies with part 15 of the FCC rules and with RSS 210 of Industry Canada. Operation is subject to the following two conditions: (1) This device does not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

<b>Guidance and manufacturer's declaration - electromagnetic emissions</b>		
The Model Smart Control is intended for use in the electromagnetic environment specified below. The customer or the user of the Model Smart Control should assure that it is used in such an environment.		
Emission test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 2	The Model Smart Control must emit electromagnetic energy in order to perform its intended function. Nearby electronic equipment may be affected.
RF emissions CISPR 11	Class B	The Model Smart Control is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable	

Table 201

## Guidance and manufacturer's declaration – electromagnetic immunity

The Model Smart Control is intended for use in the electromagnetic environment specified below. The customer or the user of the Model Smart Control should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD)  IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst  IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge  IEC 61000-4-5	±1 kV line(s) to line(s) ±2 kV line(s) to earth	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines  IEC 61000-4-11	<5% $U_T$ (>95% dip in $U_T$ ) for 0,5 cycle  40% $U_T$ (60% dip in $U_T$ ) for 5 cycles  70% $U_T$ (30% dip in $U_T$ ) for 25 cycle  <5% $U_T$ (>95% dip in $U_T$ ) for 5 sec	Not applicable	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Model Smart Control requires continued operation during power mains interruptions, it is recommended that the Model Smart Control be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field  IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE:  $U_T$  is the a.c. mains voltage prior to application of the test level.

## Guidance and manufacturer's declaration – electromagnetic immunity

The Model Smart Control is intended for use in the electromagnetic environment specified below. The customer or the user of the Model Smart Control should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
			Portable and mobile RF communications equipment should be used no closer to any part of Smart Control, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. <b>Recommended separation distance</b>
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	Not applicable	$d = 1,2 \sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m	$d = 1,2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2,3 \sqrt{P}$ 800 MHz to 2,5 GHz
			where $P$ is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and $d$ is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, <sup>a</sup> should be less than the compliance level in each frequency range <sup>b</sup> . Interference may occur in the vicinity of equipment marked with the following symbol: 

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

<sup>a</sup>Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Model Smart Control is used exceeds the applicable RF compliance level above, the Model Smart Control should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Model Smart Control.

<sup>b</sup>Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

## Recommended separation distances between portable and mobile RF communications equipment and the Model Smart Control

The Model Smart Control is intended for use in an electromagnetic environment in which radiated RF Disturbances are controlled. The customer or the user of the Model Smart Control can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Model Smart Control as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation Distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1,2 \sqrt{P}$	80 MHz to 800 MHz $d = 1,2 \sqrt{P}$	800 MHz to 2,5 GHz $d = 2,3 \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance  $d$  in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where  $P$  is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Tableau 206

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20 Beasley Drive, P.O. Box 9017, Kitchener, ON N2G 4X1 Canada

### Australie

Level 2, Norwest Quay,  
21 Solent Circuit,  
Baulkham Hills, NSW,  
NSW 2153

### Belgique

Baron de Vironlaan, 60  
b-1700 Dilbeek

### Canada

20 Beasley Drive, P.O.  
Box 9017, Kitchener, ON  
N2G 4X1

### Chine

No. 78, Qi Ming Road,  
Export Processing Zone  
District B, Suzhou  
Industrial Park, P.R.  
215021

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5 bis, rue Maryse  
Bastie - BP 15 - 69671  
Bron Cedex

### Inde

B-308, Everest Grande,  
Mahakali Caves Road,  
Andheri East,  
Mumbai-400 093

### Pays-Bas

Ir. D.S. Tuijnmanweg  
10, 4131 PN VIANEN

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10/215 Rosedale Rd.,  
M277 Private Bag,  
300987, Albany,  
Auckland

### Russie

Tverskaya str. 12,  
bldg. 9, 1th floor,  
office 98  
125009 Moscow

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First Floor –  
Selborne House,  
Fourways Golf Park,  
Roos Street,  
Fourways,  
Johannesburg

### Espagne

Avenida de la  
Industria nº13/15,  
Polígono de  
Canastell, 03690  
Sant Vicent del  
Raspeig. Alicante,  
Espagne

### Suède

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### Royaume-Uni

St. George House,  
Cygnet Court,  
Centre Park,  
Warrington,  
Cheshire WA1 1PD

### États-Unis

14755  
27<sup>th</sup> Avenue North  
Plymouth, MN 55447

Unitron Hearing Ltd.

20 Beasley Drive, P.O. Box 9017,  
Kitchener, ON N2G 4X1 Canada

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